instruction for use
gebruiksaanwijzing
mode d'emploi
gebrauchsanweisung

WE CARE U SPORT.

Fittner TOUCH & GO

percussion massage device for warming-up,
recovery and muscle relaxation

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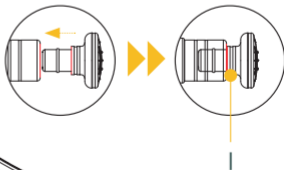
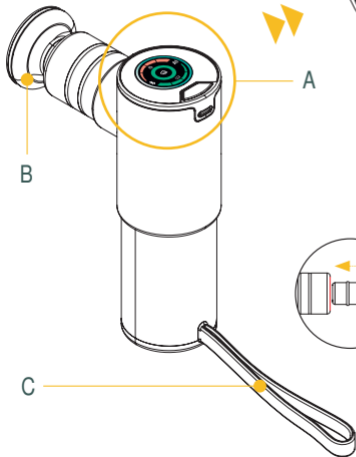
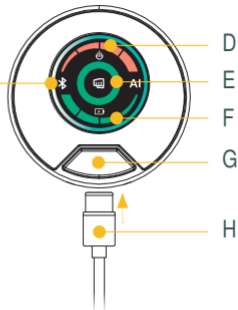
Incorrect application can cause injuries. Do not use the device if you have any physical health concerns. Read the instruction for use before using this product. Dispose of contents/container in accordance with local/regional/national/international regulations.

- | | |
|--|--------------------------------|
| A Display panel | E Pressure indicator |
| B Quick detachable massage head | F Power indicator |
| C Non-slip wrist strap | G Multi-function button |
| D Gear indicator | H Charging port |
-

1. CONTENTS

- 1x Instruction for use
 - 1x Portable Massage gun
 - 1x 2500 mAh Li-ion
 - 1x USB Type-A to Type-C charging cable
 - 4x Detachable massage attachment
-

2. DESCRIPTION

This is a Portable Massage Device with 5 powerful speed levels. It delivers high penetration muscle relief and a faster, more effective recovery for sport athletes.

Main functions:

Relieving muscle soreness, stimulating the circulation of blood, relieving stiffness of cervical vertebrae, and improving muscle recovery.

FEATURES

- 5 speed level frequency
- 2500 mAh Lithium Battery
- 1 year warranty

3. INSTRUCTIONS FOR USE

1. Press and hold the multi-function button (G) for 1 second to turn on the power.
2. Short press the multi-function button (G), the machine starts to work, each press will increase the vibration strength by 1 gear. The fifth gear is AI mode.
3. In the process of use, the pressure indicator (E) will display the current pressure in real time. There are three force gears.
4. Press and hold the multi-function button (G) for 1 second to shut down.

CHARGING

The charging port is on the top of the massage gun.

1. Please fully charge the battery for the first use.
2. Turn it off when charging.
3. Connect the power supply (H).
The battery icon will flash (F) and it shows that the machine is charging.
4. When four indicator lights (F) are on, the battery is full.
5. Do not move the machine when it is charging

REPLACE THE MESSAGE HEAD

With the machine turned off, slowly insert the massage head (B) into the barrel. As shown in the figure, the red line position of the massage head can be aligned with the red line position (I) of the machine.

Spinal head for spine massage. Strong intensity.

Round head whole body muscle massage and relaxation.

Warmy head for warm up muscle and fascia relaxation.

Bullet head suitable trigger point for massage. **DO NOT** use longer than 15 seconds on a trigger point.

MAINTENANCE

- When you do not use the product for a long time, turn it off and hold electricity while safe-keeping it. Store it in a safe, cool, and dry place.
- Charge the product every half a year if you do not use it for a long time. Holding electricity prevents battery damage due to stains.
- Please choose neutral detergent and soft cloth for daily maintenance. Do not use engine oil, organic solvents or other chemicals for cleaning.
- Please turn off the power after using the product, or before cleaning the product.
- The product should be well fixed during transportation. Avoid violent vibration, collision, rolling and falling

4. WARNINGS & PRECAUTIONS

Keep out of reach of children. Incorrect application can cause injuries. Consult a qualified person for the correct application of this device. If you are receiving medical treatment/or have unexplained pain; please ask your doctor before using this device.

DO NOT USE

If you have any physical health concerns, or if you have any of the following conditions: pacemaker usage, pregnancy, diabetes, recent surgery, epilepsy or migraine, herniated disks, spondylolisthesis or spondylosis, recent joint replacements, and metal pins/plates.

Ongoing research on the effect of vibration training toward specific medical diseases may reduce the number of contraindications. Practical experience shows that, in many cases, vibration-training-combined treatment is also an effective medical guidance. However, it must be done with doctor advice

The product is not an alternative or complement of medical treatment. If you are feeling pain, undergoing medical treatment, or using any other medical devices, please consult your doctor before using this product

CAUTION

- do not put the device into water or other liquids
- do not use this device for more than 30 minutes at a time
- do not apply this device on you head, neck or near the genitals
- do not use this device to supplement or replace any medical treatments
- do not use if you have any unexplained pain
- do not use if you are using any medical equipment

Do not touch this part when the machine is running without the massage heads being installed.

Do not touch this part when the machine is running with the massage heads installed.

5. 1-YEAR WARRANTY

If within one year from the date of purchase, this product fails due to a defect in material or workmanship, the inter-national distributor, dealer or retailer from whom the purchase was made, will repair or replace the product, or necessary components, free of charge. You must submit the proof of purchase, including the unique identification number on the box, otherwise you will be excluded from the warranty.

WARRANTY EXCLUDES

- a) damage caused by accident, abuse, mishandling or transport;
- b) units subjected to unauthorized repair (repairs carried out by persons other than authorized repairers of the producer/distributor;
- c) units not used in accordance with the device care instructions;
- d) damage exceeding the cost of the product;
- e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises;
- f) failure to provide the dated proof of purchase including the unique identification number on the box.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, and other rights may vary from country to country, and from state to state. This warranty service is available through the opening policies and procedures at Medicc Sports Europe B.V.

If you are a client, please contact the appropriate distributor, business entity, dealer or retailer from whom you directly purchased the product for warranty and follow the authorization procedures. This warranty is initiated and executed by Medicc Sports Europe B.V.

6. TECHNICAL SPECIFICATIONS

Product name: Fittner	Charging time: About 3 hrs
Type: Touch & Go	Stroke: 6 mm
Rated input: Type-c, 5 V --- 2 A	Maximum power: 70 W
Rated voltage: 12 V	Net weight: 500 g
Rated speed: 900-3200 rpm	Size: 142 x 90 x 44 mm
Battery capacity: 2500 mAh	Service temperature: 0-40 °C

7. COMPLAINTS & CUSTOMER SERVICE SUPPORT

Contact the distributor with questions/complaints or to report an adverse reaction.

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