WE CARE U SPORT.

instruction for use
gebruiksaanwijzing
mode d'emploi
gebrauchsanweisung

Fittner^{PRO}

percussion massage device for warming-up,
recovery and muscle relaxation

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

FR

DE

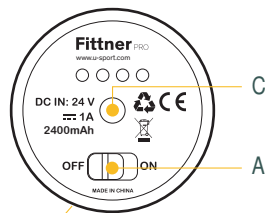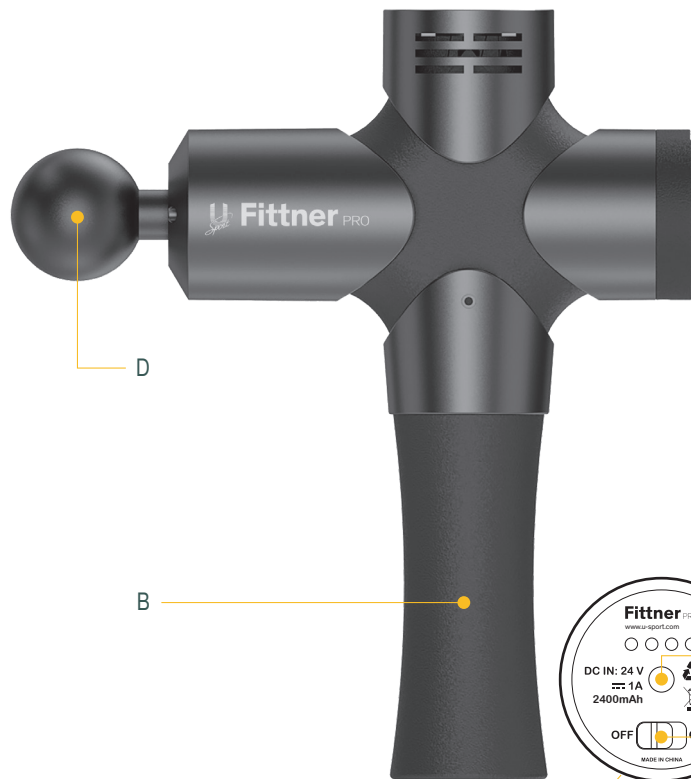Incorrect application can cause injuries. Consult a qualified person for the correct application. Please contact your physician, before using the device. Do not use the device if you have any physical health concerns. Read the instruction for use before using this product. Dispose of contents/container in accordance with local/regional/national/international regulations.

- | | |
|----------------------------------|---|
| A Power switch | E Power indicator light |
| B 2400 mAh Li-ion | F Level indicator light |
| C Charging interface | G Gear switch (speed setting button) |
| D Detachable massage head | |

1. CONTENTS

- 1x Instruction for use
- 1x U-Sport Fitter PRO massage device
- 1x 2400 mAh Li-ion
- 1x Charger (EU)
- 6x Detachable massage attachment

2. DESCRIPTION

This is a Portable Massage Device with 4 powerful speed levels. It delivers high penetration muscle relief and a faster, more effective recovery for sport athletes.

FEATURES

- noise reduction function
- 4 speed level frequency
- 2400 mAh Lithium Battery (3-6 hours battery life)
- 1 year warranty

3. INSTRUCTIONS FOR USE

Do not use this device for more than 30 minutes at a time.

This device should only be used with the provided, certified, charger.

CHARGING


Before connecting the device to the power supply, please ensure that the supply-voltage stated on the rating plate is compatible with the unit. Charge the battery before the first use.

1. It takes 12 hours for the first charge. The charging time for each subsequent time is 2 hours.
2. Make sure the device is off.
3. Connect the charger.
4. Power indicator lights, light up and starts flashing, see image 1.
5. Charging finishes when the four indicator lights shine without flashing, see image 2.

ATTACHING THE MASSAGE HEAD

1. Turn off the machine.
2. Apply the attachment into the front of the device and press firmly, see image 3.
3. When you want to change the massage head then turn and pull the head from the device, see image 4.

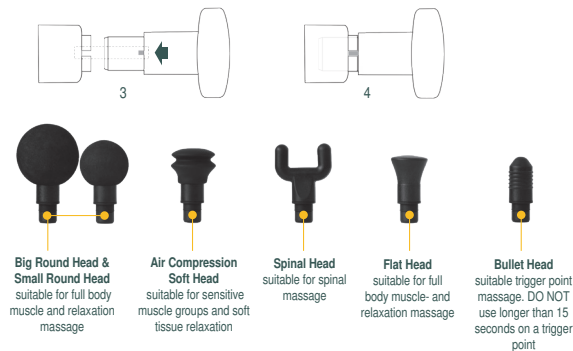START/ADJUST SPEED LEVEL

1. Turn on the power switch (A) - power indicator light will be on (E).
2. Press the Gear Switch (G) - speed setting unit will vibrate at the first gear speed.
3. To adjust the speed, press the Gear Switch.
4. The device changes speed each time you press.
5. The device automatically turns off every 10 minutes during use.
6. Press the Gear Switch again to continue using the device.
7. TURN OFF: Press the Gear Switch and turn off the Power Switch.

DO NOT

Adjust the speed / frequency (or pull out the battery) while the device is on.

CLEANING

- use a towel or a soft cloth to wipe the surface moisture
- keep the product away from heat, sources & liquids
- protect from sunlight
- dry in the shade

MAINTENANCE

- batteries in packs are not designated for exchange or removal by the users
- if problems with the power supply should occur, please contact the service center
- if the device malfunctions, please do not try to repair or to fix it by taking the device apart yourself (this leads automatically to the automatic expiry of warranty)

In case you are taking the device with you on a flight it is recommended to bring the battery in your hand-luggage.

4. WARNINGS & PRECAUTIONS

Read the instruction for use before using the product. Keep out of reach of children. Incorrect application can cause injuries. Consult a qualified person for the correct application of this device. If you are receiving medical treatment / unexplained pain; please ask your doctor before using this device.

DO NOT USE

If you have any physical health concerns or if you have any of the following conditions: pregnancy, diabetes, recent surgery, epilepsy or migraine, herniated disks, spondylolisthesis or spondylosis, recent joint replacements or IUD's and metal pins/plates

CAUTION

- do not put the device into water or other liquids
- do not use this device for more than 30 minutes at a time
- do not apply this device on you head or near the genitals
- do not use this device to supplement or replace any medical treatments
- do not use if you have any unexplained pain
- do not use if you are using any medical equipment

5. 1-YEAR WARRANTY

If within one year from the date of purchase, this product fails due to a defect in material or workmanship, the international distributor, dealer or retailer from whom the purchase was made, will repair or replace the product, or necessary components, free of charge. You must submit the proof of purchase, including the unique identification number on the box, otherwise you will be excluded from the warranty.

WARRANTY EXCLUDES

- a) damage caused by accident, abuse, mishandling or transport;
- b) units subjected to unauthorised repair (repairs carried out by persons other than authorised repairers of the producer/distributor;

- c) units not used in accordance with the device care instructions;
- d) damage exceeding the cost of the product;
- e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises;
- f) failure to provide the dated proof of purchase including the unique identification number on the box.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you. This warranty guarantees specific legal rights, and other rights may vary from country to country, and from state to state. This warranty service is available through the opening policies and procedures at Medicc Sports Europe B.V. If you are a client, please contact the appropriate distributor, business entity, dealer or retailer from whom you directly purchased the product for warranty and follow the authorization procedures. This warranty is initiated and executed by Medicc Sports Europe B.V.

TECHNICAL SPECIFICATIONS

AC Charger	
Rated input	100/240V - 50/60Hz, 24W
Rated output	24V/1A
Battery	
Rated voltage	21.6V 6INR19/66
Type	Li-ion
Capacity	2400 mAh
Working hour	3 - 6h
Charging time	120 min
Massager	
No load speed	1300-3300 min / max
Power	126W max
Frequency	1st gear: 1300 2nd gear: 1800 3rd gear: 2500 4th gear: 3300
Weight	1200g/2400g (incl. charger & retailbox)
Dimension Packaging	39 x 23,5 x 9 cm
Acoustic pressure level IpA	72dB (A)
Acoustic power level IpA	73dB (A)

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